

Softform Premier



Pressure reducing static mattress



The *Invacare* **Softform Premier** is a high specification anti-decubitus mattress designed to meet the demands of the modern ward environment. Practical and durable, this mattress delivers exceptional levels of comfort and pressure reduction.


Pressure ulcer prevention

The *Invacare* **Softform Premier** features a unique contoured insert pad which provides excellent patient comfort and optimises pressure reduction. The mattress cover is two-way stretch and waterproof, and positively helps to reduce shear and friction forces.

Cover is resistant to fluid ingress

Cover is two-way stretch waterproof polyurethane cover with high frequency welded seams; and a full length flap conceals the zip. These features combine to effectively reduce the potential of fluid ingress. Also features a toughed PU coated base to prolong the longevity of the mattress.

Features and Options

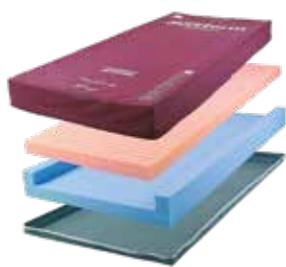

Single piece foam "U" core

Facilitates patient transfer by providing effective side-wall support.




Full length flap conceals zip

The full length flap protects the zip. This reduces the possible risk of fluid ingress, preventing contamination of the foam core.



Does not require turning or rotating

This mattress does not require turning or rotating, significantly minimizing manual handling issues for clinical personnel.

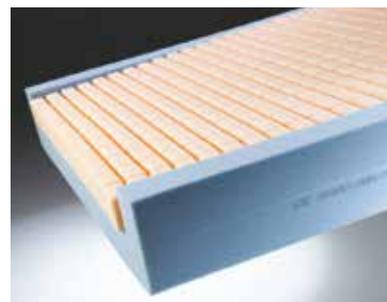

Tapered surface in the heel zone

The top surface of the mattress tapers down at the ends to relieve pressure on the heels, which are particularly vulnerable to pressure ulcers.

Technical data

For more comprehensive pre-sales information about this product, including the product's user manual, please see www.invacare-eu-export.com/softform-premier

	Width	Length	Height	Total product weight	Max. user weight	Max. cleaning temperature	Manufactured to European Standards: BS EN 597-1 Cig, BS EN 597-2 Match. Foam Interior: BS 5852: 1990 Crib Source 5; Cover BS7175: 1989 Crib Source 5.
Softform Premier	 Single size: 660 - 900 mm Double size: 910 - 1370 mm King size: 1380 - 1525 mm	 1410 - 2100 mm	 152 mm	 13 kg	 247.6 kg	 80°	

Colours

